

COMPACT MASTER ALARM

FEATURES

- 5.7" Liquid Crystal Display (LCD)
- Supply Power: 100V-240V AC and Circuit breaker protected
- Alarm Point Inputs: up to 30 points
- Remote alarm communications: Modbus (standard); BACnet (optional); Ethernet (optional)
- NFPA[®] 99 Compliant
- Manufactured in an ISO 9001 and 13485 Environment
- Enclosure materials constructed of heavy duty steel to protect internal components from adverse conditions and electrical interference
- Alarm conditions identified via audible alarms and visual indicators
- Audible Alarm-Point Silence Button provided to silence the alarm (field adjustable duration)
- Audible Alarm Horn is adjustable
- Audible Alarm is to remain in an alarm state until the condition that initiated the alarm condition is corrected
- Visual Alarm-Point Indicators shall illuminate and/or flash even if the Audible Alarm is silenced and, to remain in an alarm state until the condition initiating the alarm is corrected




PANEL SETTINGS

- Completely field adjustable and password protected
- All settings accomplished via front screen panel
- Selectable English or Spanish Instructions
- Password Protected (Customizable)
- Field Adjustable 0-30pt alarm set-up
- Adjustable N.O./N.C. input selection with Set-All feature
- All Alarm Points in an alarm condition are viewed in priority scrolling
- Pre-configured or customizable Alarm Point names
- Auto scroll with manual scroll options for non-alarm condition viewing
- Maintenance Mode for set-up and alarm point maintenance



PRODUCT DIMENSIONS

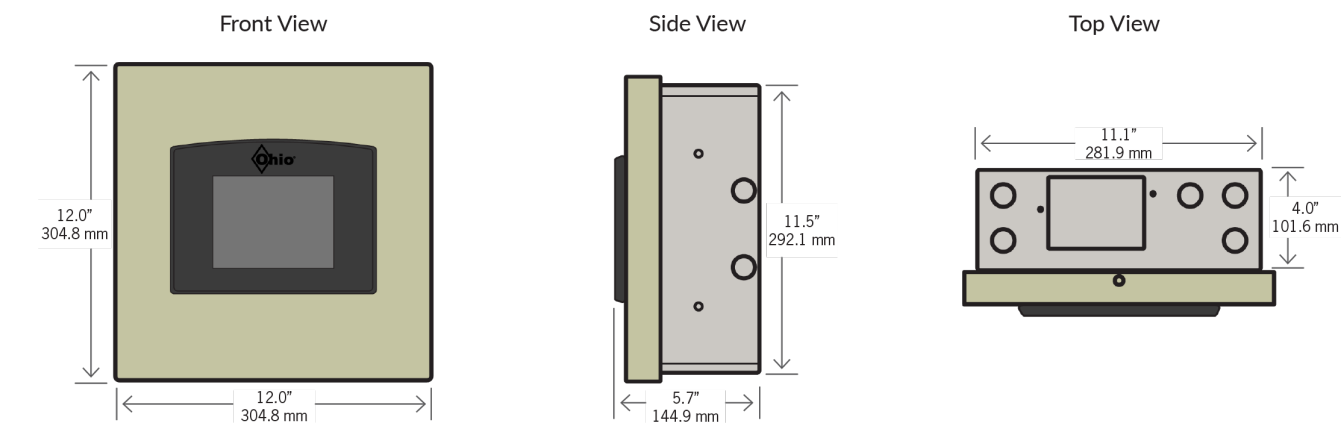

ORDERING INFORMATION

Available with Modbus, BACnet or Ethernet. All units come with a 5.7" full color LCD touchscreen and can be set up for one to thirty unique alarm points inputs

Part Number	Description*
263859	Compact Master Alarm HMI 30 Point
263859-B	Compact Master Alarm HMI 30 Point with BACnet
263859-E	Compact Master Alarm HMI 30 Point with Ethernet


*Modbus included with every unit.

PRESSURE/VACUUM PIPEAWAY ASSEMBLIES WITH SWITCHES

Low Pressure Switch

Part Numbers


- 263741 O₂
- 263743 Air
- 263744 N₂O
- 263747 CO₂
- 263749 INSTR AIR - L

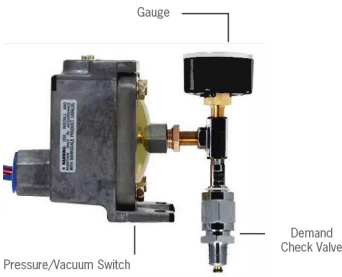

High Pressure Switch

Part Numbers

- 263746 N₂
- 263746 -IA INSTR AIR - H

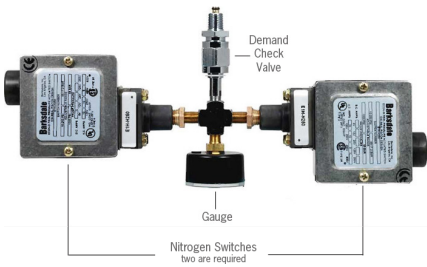

Vacuum Pressure Switch

Part Numbers

- 263742 Vac
- 263745 WAGD



www.gce-medical.com
www.ohiomedical.com